
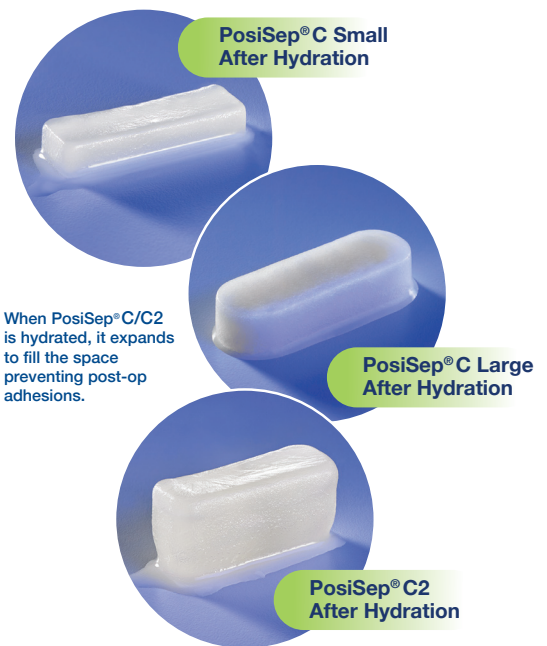


For Hospital, ASC, and Office Based Procedures

PosiSep[®] C/C² are engineered for use in Hospital, ASC, and Office Based procedures, and are intended to minimize bleeding and edema, prevent adhesions between the septum and nasal cavity, and act as an adjunct to aid in the natural healing process. PosiSep[®] C/C² offer a slender, easily placed, nasal dressing which allow gentle placement and rapid expansion conforming to the anatomy upon hydration. PosiSep[®] C/C² are constructed of a patient-comfortable sponge which is manufactured from naturally occurring chitosan polymers. PosiSep[®] C/C² are available in three sizes to choose from for various surgical applications. PosiSep[®] C² expands to a greater volume (0.75") to fill larger spaces of a sinus cavity.

Designed for patient comfort and simple application.

- ✓ Widely studied and published anti-adhesion properties of our proprietary NOCC technology
- ✓ Biocompatible – made from naturally occurring chitosan polymers
- ✓ Aids in the natural wound healing process
- ✓ Easy to apply-no mixing or additional prep time necessary
- ✓ Eliminated within 4-5 days by way of daily irrigation
- ✓ No special storage required
- ✓ Rapidly expands upon hydration
- ✓ Offers easier placement with superior field of view

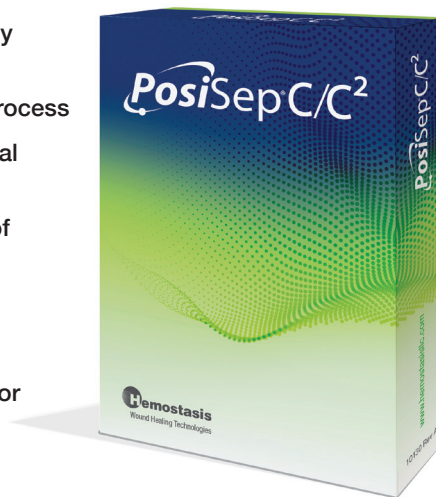

ORDERING INFORMATION: PosiSep[®] C/C² is available in three (3) sizes. See ordering information below.

Part Number:	Description:	How Supplied:
9210754	PosiSep [®] C Nasal Dressing 0.5" x 1.5" 5pk (US)	PosiSep [®] C & C ² are packaged five (5) sterile, single use units per carton. One (1) sterile unite contains one (1) device. PosiSep [®] C & C ² are terminally sterilized by gamma irradiation.
9210854	PosiSep [®] C Nasal Dressing 0.6" x 2" 5pk (US)	
9210764	PosiSep [®] C ² Nasal Dressing 0.5" x 1.5" 5pk (US)	

CAUTION: Federal (U.S.) law restricts this device to sale by or on the order of a physician. PosiSep[®] C/C² is contraindicated for patients who have known allergies to shellfish. For detailed information regarding indications for use, warnings and precautions, see Instructions for Use.

To place an order or for more information on PosiSep[®] C/C² or our other medical technologies, please contact a sales representative by calling our toll free number 866-612-2568