



Designed specifically for providing localized drug delivery in the sinus frontal ostia.

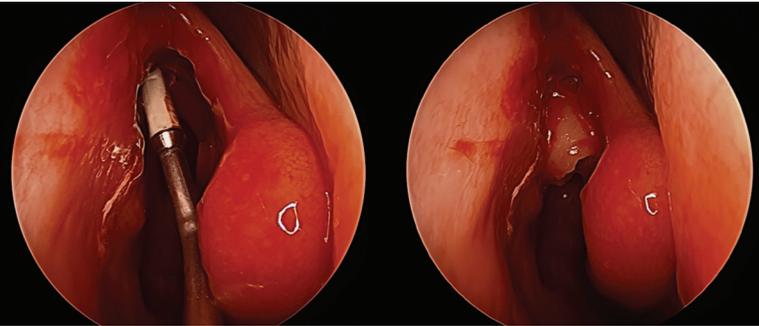

VECTRA™ puts localized drug delivery in the hands of surgeons.

VECTRA™ was specifically developed to administer medicinal substances to treat nasal disorders. The patented VECTRA™ microsponge is supplied in a compressed tubular design preloaded on one end of an insertion tool. The opposite end of the tool has a Luer fitting which is connected to a syringe containing the medicine to be delivered. The tool is used to place the VECTRA™ microsponge in the frontal ostia. Once in place, the medicine can be delivered to the VECTRA™ device. As the medicine hydrates the VECTRA™ microsponge, the microsponge expands to conform to the anatomy of the ostia. The tool is then removed and the VECTRA™ microsponge remains in place. The intimate contact allows for localized drug delivery to the desired target area.

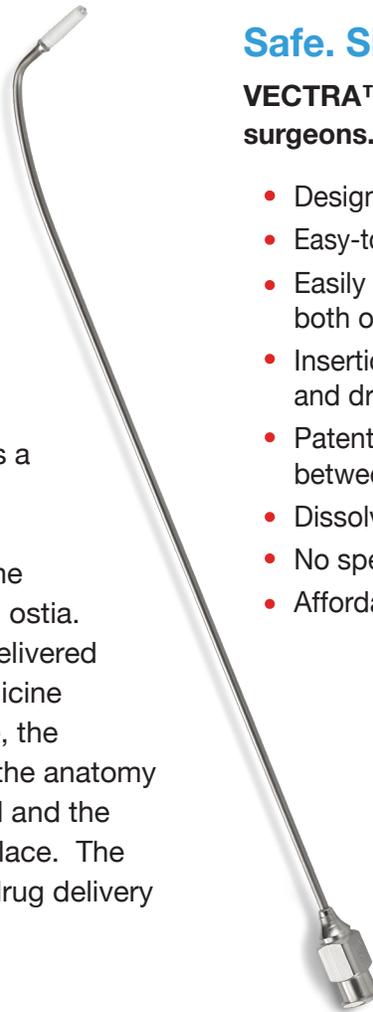

Safe. Simple. Effective.

VECTRA™. Designed in conjunction with sinus surgeons.

- Designed for localized drug delivery
- Easy-to-use
- Easily placed and hydrated in the frontal ostia in both operating room and clinic
- Insertion tool provides easy, accurate placement; and drug delivery to microsponge
- Patented tubular design maintains an open lumen between the sinus and nasal cavity
- Dissolves away in 3-6 days
- No special storage or refrigeration required
- Affordable

ORDERING INFORMATION

| Part Number: | Description: | How Supplied: |
|--------------|-------------------------------------|---|
| 9210864 | VECTRA™ F drug delivery microsponge | VECTRA™ F drug delivery device is packaged with two (2) sterile, single use units per box. One (1) sterile unit contains one (1) VECTRA™ microsponge and one (1) frontal sinus insertion tool. VECTRA™ is terminally sterilized by gamma irradiation. |

CAUTION: Federal (U.S) law restricts this device to sale by or on the order of a physician. For detailed information regarding indications for use, warnings and precautions, see Instructions for Use.

Kenalog is a registered trademark of Westwood-Squibb Pharmaceuticals Inc.

To place an order or for more information on VECTRA™ F or our other medical technologies, please contact a sales representative by calling our toll free number 866-612-2568